




Simplicity when you want it. Security when you need it.

Designed to provide a safe, efficient, and reproducible solution for anterior cervical plating

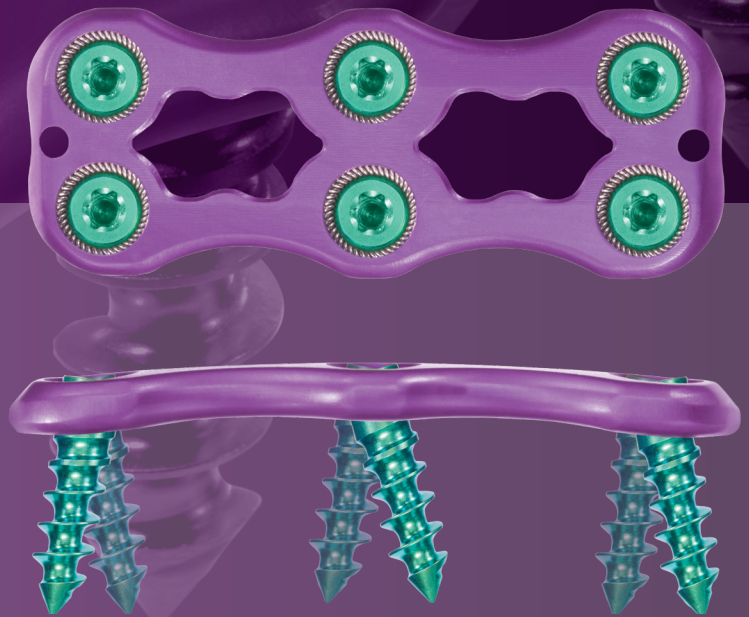

NuVasive Helix ACP Overview


NuVasive Helix ACP

- NuVasive Helix ACP plate has a 2.5mm profile to reduce the chance of postoperative discomfort.
- 5-Level plates available via special order.
- Fixed- and Variable-Angle Screws for the perfect construct type, depending on patient needs.
- Variable-Angle Screws feature a 20° cone-of-angulation to accommodate surgeon needs.



NuVasive Helix ACP Canted Coil Lock (CCL)

- Consistently and securely blocks screws into the plate while still allowing the plate to be securely lagged to bone for a tight fit.
- Screws can be removed and repositioned easily while maintaining CCL integrity.
- Aggressive screws in a full complement of lengths and diameters accommodate a variety of surgical scenarios.

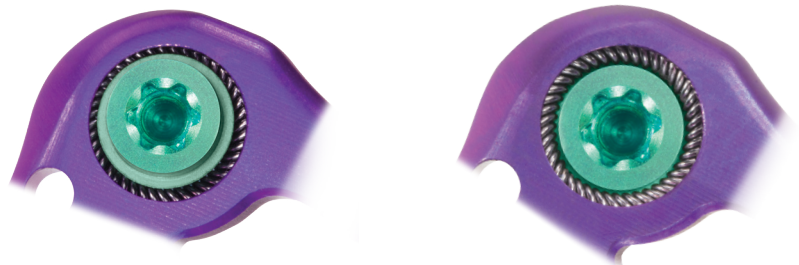

UNLOCKED


Ledge of screw is visible, and head of screw is proud to plate surface.

LOCKED

Ledge of screw is completely covered, and head of screw is seated below NuVasive Helix ACP Canted Coil Lock.

Balanced and Intuitive Instrumentation


- DTS Guides allow pre-drilling, tapping, and screw placement without changing insertion angle.

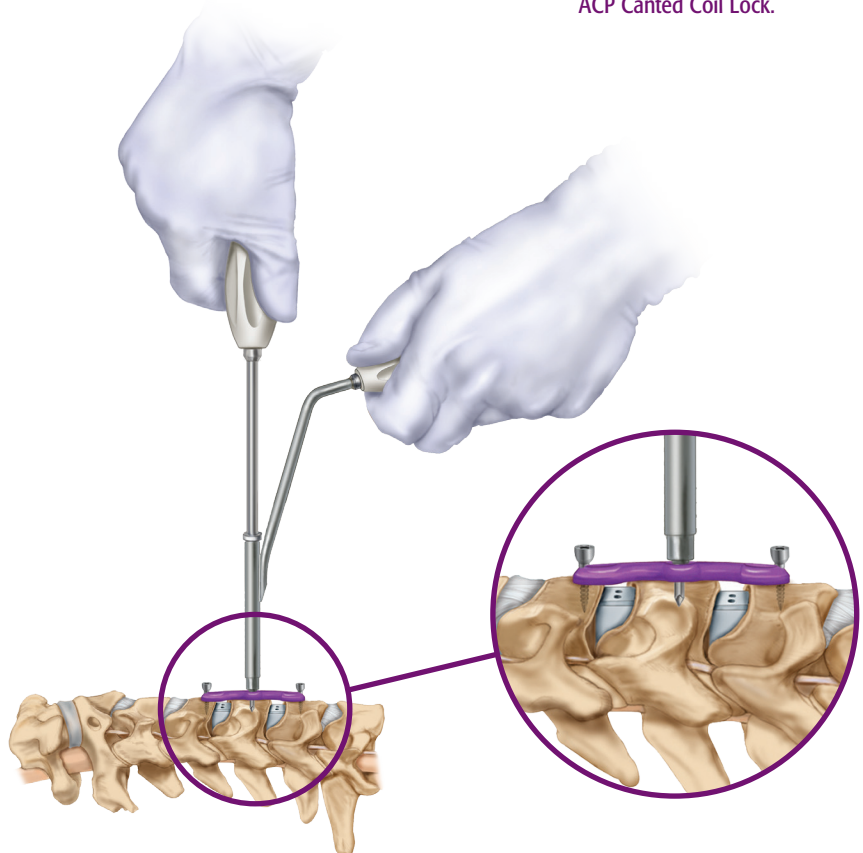

NuVasive Helix ACP Overview

Low-Profile Leading Edge

- Plate sits flush against anterior cervical spine.
- Decreased “step” on leading edge of plate to reduce esophageal discomfort postoperatively.
- 1mm leading edge.

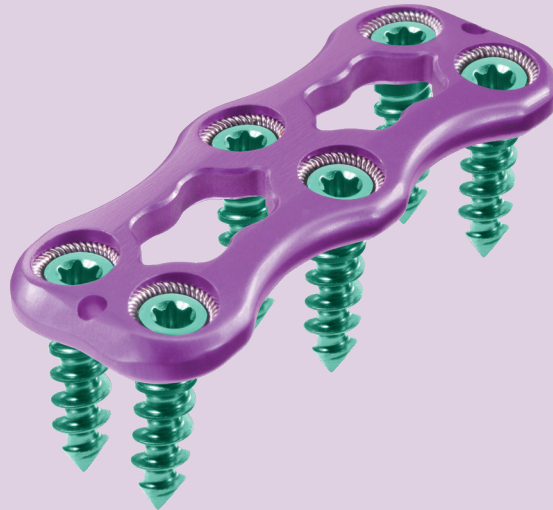

Narrow Plate Reduces Retraction During Placement

- Reduced retraction for easy placement.
- Less than 16mm at widest point.



Intraoperative Visualization

- Large graft windows allow for graft placement verification.



NUVASIVE HELIX ACP INSTRUMENTS

| DESCRIPTION | CATALOG # |
|---|-----------|
| 13mm Drill (Sterile) | 7730016 |
| 13mm Tap Bit | 7730026 |
| Screwdriver | 7730041 |
| Awl | 7730062 |
| Self-Centering Awl (Variable) | 7730063 |
| Self-Centering Awl (Fixed) | 7730064 |
| Screw Extractor | 7730072 |
| Drill Guide, Fixed | 7730081 |
| Drill Guide, Variable | 7730091 |
| Plate Bender | 7730701 |
| NuVasive Helix ACP Sterilization Tray | 7730002 |
| DTS Guide - Fixed, 0° | 7730086 |
| DTS Guide - Fixed, 7° | 7730088 |
| Plate Holder | 7730043 |
| Universal Handle | 7730600 |
| TEMPORARY FIXATION | |
| Temporary Fixation Pin (Temporary Tack) | 7737110 |

NUVASIVE HELIX ACP IMPLANTS

| DESCRIPTION | CATALOG # |
|---|-----------|
| SELF-DRILLING SCREWS (VARIABLE, 4.0) | |
| 4.0 x 11mm Self-Drilling Variable Screw | 7739111 |
| 4.0 x 13mm Self-Drilling Variable Screw | 7739113 |
| 4.0 x 15mm Self-Drilling Variable Screw | 7739115 |
| SELF-TAPPING SCREWS (VARIABLE, 4.0) | |
| 4.0 x 11mm Self-Tapping Variable Screw | 7733111 |
| 4.0 x 13mm Self-Tapping Variable Screw | 7739213 |
| 4.0 x 15mm Self-Tapping Variable Screw | 7739215 |
| 4.0 x 17mm Self-Tapping Variable Screw | 7739217 |
| 4.0 x 19mm Self-Tapping Variable Screw | 7733119 |
| SELF-TAPPING SCREWS (VARIABLE, 4.5) | |
| 4.5 x 13mm Self-Tapping Variable Screw | 7735113 |
| 4.5 x 15mm Self-Tapping Variable Screw | 7735115 |
| 4.5 x 17mm Self-Tapping Variable Screw | 7735117 |
| 4.5 x 19mm Self-Tapping Variable Screw | 7735119 |
| SELF-DRILLING SCREWS (FIXED, 4.0) | |
| 4.0 x 11mm Self-Drilling Fixed Screw | 7732111 |
| 4.0 x 13mm Self-Drilling Fixed Screw | 7732113 |
| 4.0 x 15mm Self-Drilling Fixed Screw | 7732115 |

NUVASIVE HELIX ACP IMPLANTS

| DESCRIPTION | CATALOG # |
|---|-----------|
| SELF-TAPPING SCREWS (FIXED, 4.0) | |
| 4.0 x 11mm Self-Tapping Fixed Screw | 7734111 |
| 4.0 x 13mm Self-Tapping Fixed Screw | 7734113 |
| 4.0 x 15mm Self-Tapping Fixed Screw | 7734115 |
| 4.0 x 17mm Self-Tapping Fixed Screw | 7734117 |
| 4.0 x 19mm Self-Tapping Fixed Screw | 7734119 |
| SELF-TAPPING SCREWS (FIXED, 4.5) | |
| 4.5 x 13mm Self-Tapping Fixed Screw | 7736113 |
| 4.5 x 15mm Self-Tapping Fixed Screw | 7736115 |
| 4.5 x 17mm Self-Tapping Fixed Screw | 7736117 |
| PLATES | |
| 1-Level NuVasive Helix ACP, 20mm | 7738120 |
| 1-Level NuVasive Helix ACP, 22mm | 7738122 |
| 1-Level NuVasive Helix ACP, 24mm | 7738124 |
| 1-Level NuVasive Helix ACP, 26mm | 7738126 |
| 1-Level NuVasive Helix ACP, 28mm | 7738128 |
| 1-Level NuVasive Helix ACP, 30mm | 7738130 |
| 1-Level NuVasive Helix ACP, 32mm | 7738132 |
| 1-Level NuVasive Helix ACP, 34mm | 7738134 |
| 2-Level NuVasive Helix ACP, 36mm | 7738236 |
| 2-Level NuVasive Helix ACP, 38mm | 7738238 |
| 2-Level NuVasive Helix ACP, 40mm | 7738240 |
| 2-Level NuVasive Helix ACP, 42mm | 7738242 |
| 2-Level NuVasive Helix ACP, 44mm | 7738244 |
| 2-Level NuVasive Helix ACP, 46mm | 7738246 |
| 2-Level NuVasive Helix ACP, 48mm | 7738248 |
| 2-Level NuVasive Helix ACP, 50mm | 7738250 |
| 2-Level NuVasive Helix ACP, 52mm | 7738252 |
| 2-Level NuVasive Helix ACP, 54mm | 7738254 |
| 3-Level NuVasive Helix ACP, 56mm | 7738356 |
| 3-Level NuVasive Helix ACP, 58mm | 7738358 |
| 3-Level NuVasive Helix ACP, 60mm | 7738360 |
| 3-Level NuVasive Helix ACP, 62mm | 7738362 |
| 3-Level NuVasive Helix ACP, 64mm | 7738364 |
| 3-Level NuVasive Helix ACP, 66mm | 7738366 |
| 3-Level NuVasive Helix ACP, 68mm | 7738368 |
| 4-Level NuVasive Helix ACP, 70mm | 7738470 |
| 4-Level NuVasive Helix ACP, 74mm | 7738474 |
| 4-Level NuVasive Helix ACP, 78mm | 7738478 |
| 4-Level NuVasive Helix ACP, 82mm | 7738482 |
| 4-Level NuVasive Helix ACP, 86mm | 7738486 |
| 4-Level NuVasive Helix ACP, 90mm | 7738490 |

5-Level NuVasive Helix ACP plates are not part of the standard configuration, but may be ordered separately by contacting NuVasive Customer Service.



To order, please contact your NuVasive Sales Consultant or Customer Service Representative today at:

NuVasive, Inc. 7475 Lusk Blvd., San Diego, CA 92121 • phone: 800 475 9131

NuVasive Netherland B.V. Jachthavenweg 109 A, 1081 KM Amsterdam, The Netherlands • phone: +33 20 72 33 000

nuvasive.com

©2020. NuVasive, Inc. All rights reserved. 9500171 D