



NuVasive Helix Mini ACP



The premium features you've come to expect – made faster, smaller, and more secure.

NuVasive Helix Mini ACP System Overview


NuVasive Helix Mini ACP

- Fixed- and Variable-Angle Bone Screws for the perfect construct type, depending on patient needs.
- Variable-Angle Bone Screws feature a 20° cone-of-angulation to accommodate surgeon needs.
- NuVasive Helix Mini ACP Plate has a 2.4mm profile to reduce the chance of postoperative discomfort.
- Plate sizes accommodate 1- and 2-Level construct needs.



UNLOCKED

Ledge of Bone Screw is visible, and head of Screw is proud to Plate surface.

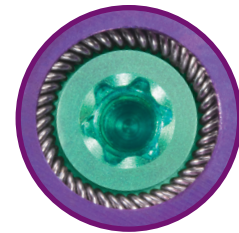

LOCKED


Ledge of Bone Screw is completely covered, and head of Screw is seated below Helix Canted Coil Lock.

Helix Canted Coil Lock (CCL)

- Helix Canted Coil Lock consistently and securely blocks Bone Screws into the Plate while still allowing the Plate to be securely lagged to bone for a tight fit.
- Bone Screws can be removed and repositioned easily while maintaining CCL integrity.
- Aggressive Bone Screws in a full complement of lengths and diameters accommodate a variety of surgical scenarios.

Balanced and Intuitive Instrumentation


- DTS Guides allow pre-drilling, tapping, and Bone Screw placement without changing insertion angle.

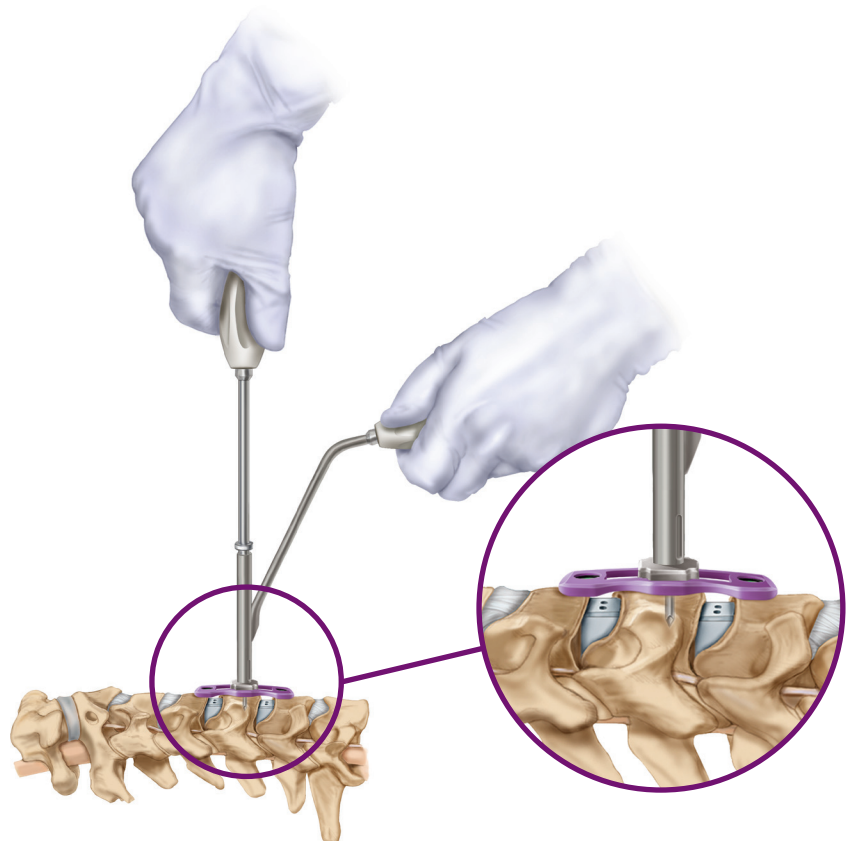

NuVasive Helix Mini ACP System Overview


Easy and Precise Intraoperative Plate Bending

- The plate may be bent intraoperatively into either further lordosis or into kyphosis in order to accommodate patient pathology.

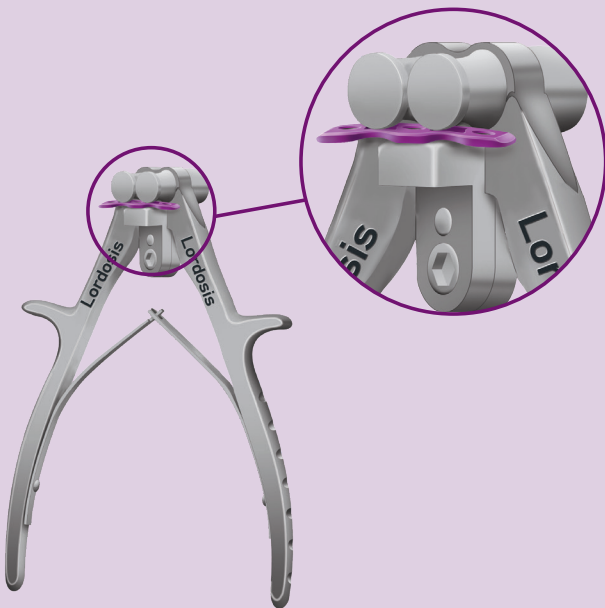


Low-Profile Leading Edge

- Plate sits flush against anterior cervical spine.
- Decreased "step" on leading edge of Plate to reduce esophageal discomfort postoperatively.

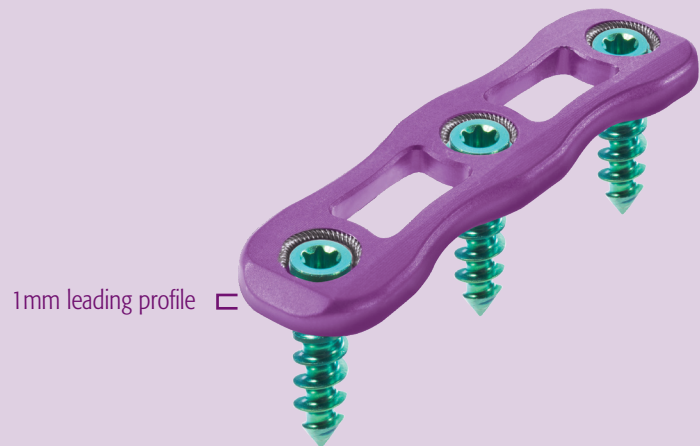

Narrow Plate Reduces Retraction During Placement

- Reduced retraction for easy placement.
- Less than 12mm at widest point.
 - 10mm at the waist.

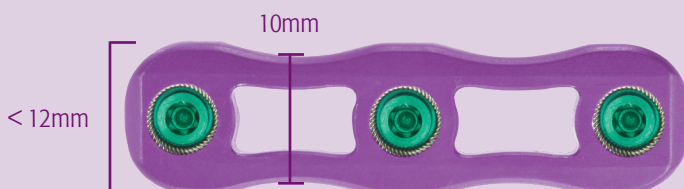

Intraoperative Visualization

- Large portals allow for graft placement verification.
- Unique graft-viewing/alignment window gives accurate reference for straight midline placement.



NUVASIVE HELIX MINI ACP IMPLANTS AND INSTRUMENTS

IMPLANTS	CATALOG #
1-Level NuVasive Helix Mini ACP, 20mm	7740120
1-Level NuVasive Helix Mini ACP, 22mm	7740122
1-Level NuVasive Helix Mini ACP, 24mm	7740124
1-Level NuVasive Helix Mini ACP, 26mm	7740126
1-Level NuVasive Helix Mini ACP, 28mm	7740128
1-Level NuVasive Helix Mini ACP, 30mm	7740130
1-Level NuVasive Helix Mini ACP, 32mm	7740132
1-Level NuVasive Helix Mini ACP, 34mm	7740134
2-Level NuVasive Helix Mini ACP, 36mm	7740236
2-Level NuVasive Helix Mini ACP, 38mm	7740238
2-Level NuVasive Helix Mini ACP, 40mm	7740240
2-Level NuVasive Helix Mini ACP, 42mm	7740242
2-Level NuVasive Helix Mini ACP, 44mm	7740244
2-Level NuVasive Helix Mini ACP, 46mm	7740246
2-Level NuVasive Helix Mini ACP, 48mm	7740248
2-Level NuVasive Helix Mini ACP, 50mm	7740250

INSTRUMENTS	CATALOG #
Bone Screw Driver	7730042
Drill Guide Awl	7730062
Self-Centering Awl (Variable)	7730063
Self-Centering Awl (Fixed)	7730064
Bone Screw Extractor	7730072
Drill Guide, Fixed Angle	7730081
Drill Guide, Variable Angle	7730091
Plate Bender	7740701
NuVasive Helix Mini ACP Sterilization Tray	7740001
7° End DTS Guide	7740086
0° Mid DTS Guide	7740088
NuVasive Helix Mini ACP Plate Holder	7740043
Universal Handle	7730600
NuVasive Helix Mini ACP - IFU	9400236
Drill Bit	7730016
Tap	7730026

IMPLANTS	CATALOG #
Self-Drilling Bone Screws (Variable, 4.0)	
4.0 x 11mm Self-Drilling Variable Screw	7741111
4.0 x 13mm Self-Drilling Variable Screw	7741113
4.0 x 15mm Self-Drilling Variable Screw	7741115
Self-Drilling Bone Screws (Fixed, 4.0)	
4.0 x 11mm Self-Drilling Fixed Screw	7742111
4.0 x 13mm Self-Drilling Fixed Screw	7742113
4.0 x 15mm Self-Drilling Fixed Screw	7742115
Self-Tapping Bone Screws (Variable, 4.0)	
4.0 x 11mm Self-Tapping Variable Screw	7748111
4.0 x 13mm Self-Tapping Variable Screw	7748113
4.0 x 15mm Self-Tapping Variable Screw	7748115
Self-Tapping Bone Screws (Fixed, 4.0)	
4.0 x 11mm Self-Tapping Fixed Screw	7747111
4.0 x 13mm Self-Tapping Fixed Screw	7747113
4.0 x 15mm Self-Tapping Fixed Screw	7747115
Self-Tapping Bone Screws (Variable, 4.5)	
4.5 x 11mm Self-Tapping Variable Screw	7735111
4.5 x 13mm Self-Tapping Variable Screw	7735113
4.5 x 15mm Self-Tapping Variable Screw	7735115
Self-Tapping Bone Screws (Fixed, 4.5)	
4.5 x 11mm Self-Tapping Fixed Screw	7746111
4.5 x 13mm Self-Tapping Fixed Screw	7746113
4.5 x 15mm Self-Tapping Fixed Screw	7746115



■ NuVasive, Inc. 7475 Lusk Blvd., San Diego, CA 92121 USA • phone: 800 475 9131

EC REP NuVasive Netherlands B.V. Jachthavenweg 109A, 1081 KM Amsterdam, The Netherlands • phone: +33 20 72 33 000

nuvasive.com

©2020. NuVasive, Inc. All rights reserved. 9500188 C

CE 2797